



Tilt Table

Head-up Tilt Testing


For the clinical evaluation of patients presenting with syncopal problems.

Designed specifically for autonomic testing, with a rapid tilt (0 to 70° in 9 seconds) and a wide choice of patient positions, including a -12° Trendelenberg.

Three pre-set tilt angles can be programmed in to the pneumatic tilt control hand switch, and then accessed at the push of a button.

- Fast and smooth transition with an electromechanical drive.
- Floor level removable footboard for stress free patient access.
- Includes 2 secure restraining belts for patient comfort and safety.
- Ergonomically designed with both clinician and patient in mind.

Model S1

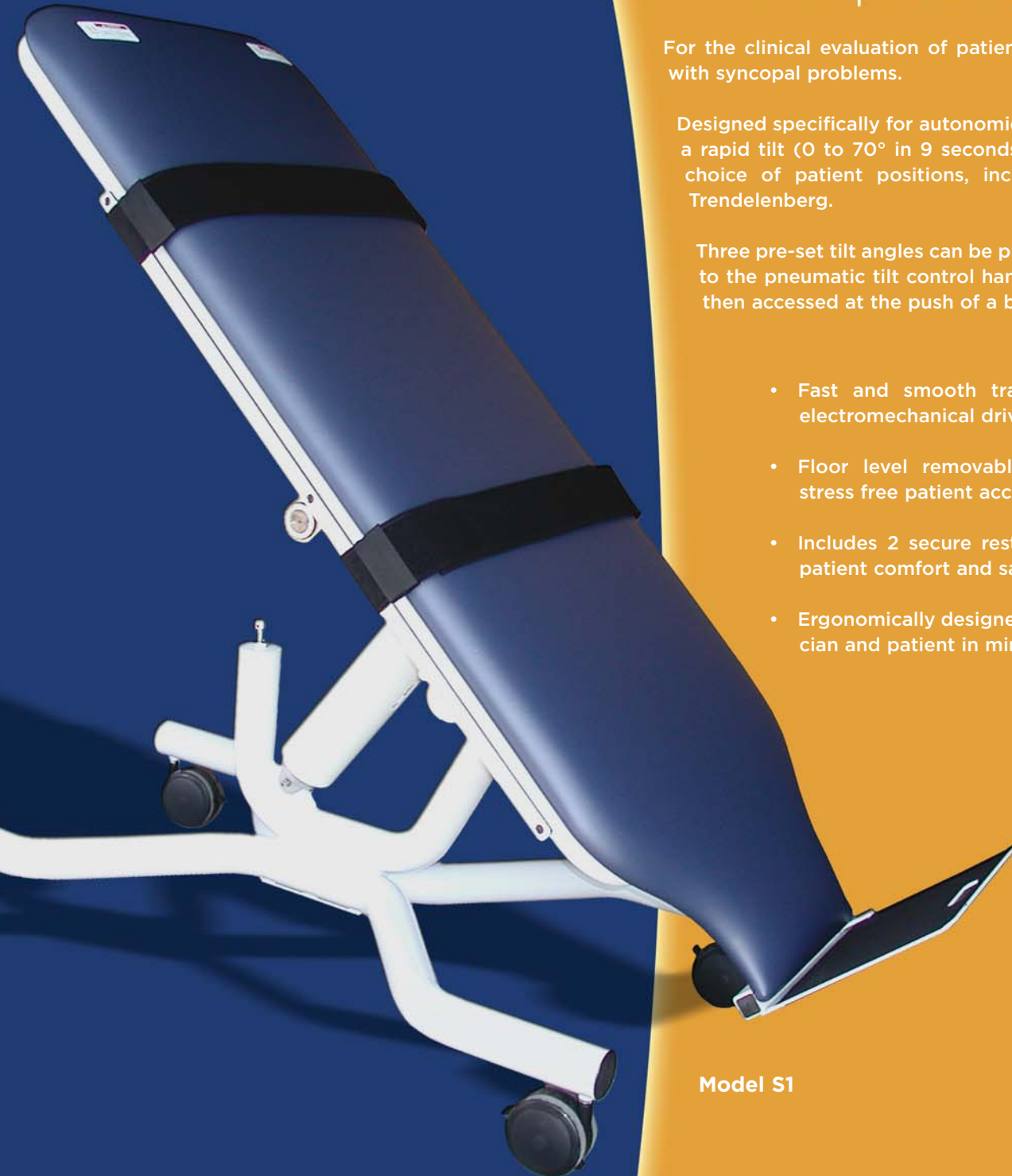

TILT TABLE TECHNICAL SPECIFICATIONS

Length	78.0"	Finish	Polyurethane Pearl Gray
Width	28.0"	Armrest	One adjustable arm assembly standard.
Height	33.0"	Electrical	TILT: Electromechanical drive, 7.9 amps, 120V, 60Hz current, equipped with hospital grade plug. Electrical components UL listed or UL recognized. (220 volt, 50Hz optional.)
Tilt Range	-12 to 85 degrees	Warranty	Two year limited warranty for structural frame and major moving parts. One year limited warranty for upholstery and minor moving parts.
Tilt Speed	9 seconds (0 to 70 degrees)	Table Load	Maximum weight capacity 350 lbs
Frame	Lower Frame (steel round tube) 3" x 11 GA and 7 GA Upper frame (steel round tube) 1 1/2" x 11 GA	Weight	Table 275 lbs Table shipping weight 350 lbs
Footboard	Removable		
Casters	Four - 5" locking		
Table Top Material	Vinyl covered 2" foam cushion over non-translucent substructure		
Table Cover Pad	Fire Retardant vinyl cover pad. Imperial Blue in color.		
Restraining Belts	Two restrainer belts with hook and loop fasteners		

Front View



At -12 degrees



At 85 degrees



With WR Testworks Autonomic Lab

For more information on this product, please contact:

WR

WR Medical Electronics Co.

www.wrmed.com

123 North 2nd Street

Stillwater, MN 55082

800.635.1312 | Fax: 800.990.9733